


Liberty™ Series – Heel Comfort

Heel Flotation Cushion

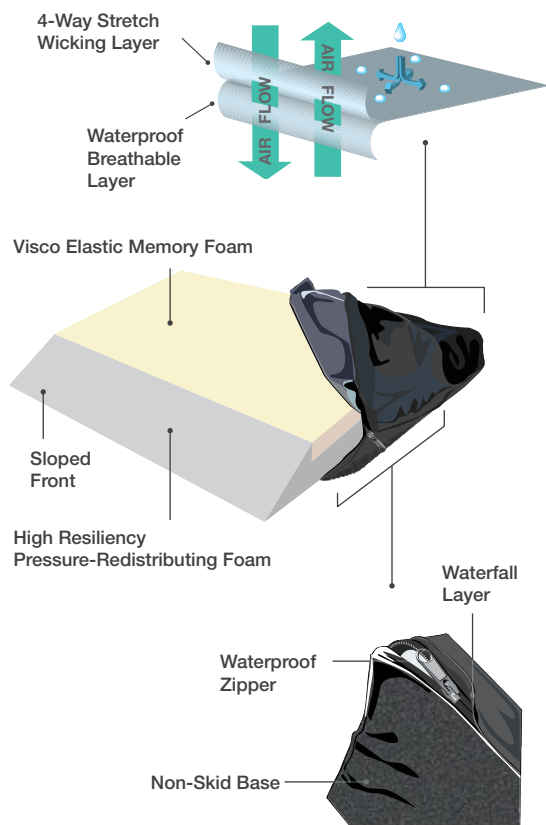

- Use positioning devices to reduce pressure from the heels and to prevent direct contact between bony prominences.
- Pillows are a temporary solution, but are often not effective because they can flatten, are not pressure relieving and they can roll off the bed or drift out of position.
- Heel flotation cushions suspend the heels while providing pressure relief by having a wider surface area for individuals who change position of his or her extremities when in bed. Make sure the cushion's cover does not stick to the skin or increase the shearing forces.
- Heel protectors with straps can be useful, but must be checked and monitored frequently to make sure straps are not too tight or loose which can lead to potential for skin breakdown. They also may cause additional heat which is one of the factors that can affect skin integrity.

Liberty™ Series – Heel Comfort

Heel Flotation Cushion

Product Specifications

Wound Risk Level: Prevention

Cushion Thickness: 3"

4"

Product Description

The Liberty Heel Comfort is designed to offer comfort and support to the lower extremities while suspending the heels. Visco elastic memory foam conforms to calves and bony prominences while floating the heels. The sloped front design supports the knees to help prevent hyper-extension. Cover is labeled with a "Float Heels This End" tag for proper use.

Models and Sizes Available

Model #	Size (WIDTH x DEPTH x HEIGHT)
H3 – 30 x 18	30" x 18" x 3"
H3 – 30 x 20	30" x 20" x 3"
H3 – 30 x 24	30" x 24" x 3"
H4 – 30 x 20	30" x 20" x 4"
H4 – 30 x 24	30" x 24" x 4"

Comfort • Safety • Mobility

Toll-Free 1(866)330-5336

Keen

Copyright © 2007, Keen Mobility Company. All rights reserved. Rev 070312. All trademarks are property of their respective owner.

Bed Safety 29